




3914 DEL AMO BLVD.
SUITE 901
TORRANCE, CA 90503
www.celprogen.com

Phone: 310 542 8822
Fax: 310 542 8028
Email: info@celprogen.com
stemcells@celprogen.com

DATA SHEET

Human Bone Marrow Derived Endothelial Stem Cells – T25 Plated Cells

Catalog number:	36096-26-T25
Description:	Human Bone Marrow Derived Endothelial Stem Cells T25 Plated Cells. Also available in T75, T150, and T225 tissue culture flask with plated Cells, shipped at room temperature or in Frozen Vial, shipped in Dry-ice. The Clonal Bone Marrow Derived Endothelial Stem Cells was derived from Adult Bone Marrow Derived Endothelial Cells. They were maintained in Celprogen's Human Bone Marrow Derived Endothelial Stem Cells Complete Growth Medium and sub-cultured every 24 to 48 hours on Human Bone Marrow Derived Endothelial Stem Cells Extra-cellular Matrix.
Source:	Adult Bone Marrow Derived Endothelial Cells
Donors:	All donors from which the Cells were derived were pre-screened; donors tested negative for the usual blood donation infectious disease panel (ABO/RH, Hepatitis B Surface Antigen, HIV1 and 2, Syphilis, hepatitis B Core, Human T Lymphocyte Virus 1 and 2, Hepatitis C Virus, Antibody Screen, Nucleic Amplification Test for HIV 1 HCV, West Nile Virus and Antibodies to Trypanosoma cruzi (the agent of Chagas disease).
Mycoplasma test:	Negative-PCR and mycoplasma agar methods
Sterility:	Negative for bacteria, yeast, and mold
Storage Conditions:	
Frozen Vial:	Liquid nitrogen vapor phase for frozen Ampule of Human Bone Marrow Derived Endothelial Stem Cells
Plated Cell Culture:	Temperature 37°C in 5% CO ₂ humidified incubator. For plated cells in tissue culture flask, upon receipt of the cells wipe the flask with 70% ethanol and transfer to sterile tissue culture hood. In the tissue culture hood remove the media from the cells and wash the cells with 1X PBS sterile solution, for 2-3 minutes, remove the 1X PBS solution and then Trypsinize. After Trypsinization of the Cells neutralize the Trypsin with equal volume of Human Bone Marrow Derived Endothelial Stem Cells Complete Growth Media with Serum and collect the Cell suspension in sterile conical centrifuge tube in the tissue culture hood. Centrifuge the cell suspension at 100g for 7 minutes in centrifuge. Plate cells 5x10 ⁵ cells per pre-coated flasks with Human Bone Marrow Derived Endothelial Stem Cells Extra-cellular Matrix for Expansion in

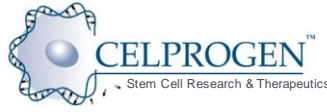

3914 DEL AMO BLVD.
SUITE 901
TORRANCE, CA 90503
www.celprogen.com

Phone: 310 542 8822
Fax: 310 542 8028
Email: info@celprogen.com
stemcells@celprogen.com

Human Bone Marrow Derived Endothelial Stem Cells Complete Growth Medium.

Positive Markers: CD15 (SSEA-1), CD90, CD105, CD106, CD117 & CD309

Morphology & Proliferation:

Mixed population of cells with approximately 95% attached cells and the other 5.0% in suspension, need to change cell culture media every day after 48 hours of initial cell culture or when the media starts changing color to slight yellow from pink. Fast growing cell culture, change media with Celprogen's Human Bone Marrow Derived Endothelial Stem Cells Complete Growth Medium with the appropriate Human Bone Marrow Derived Endothelial Stem Cells Extra-cellular Matrix. Temperature 37°C in 5% CO₂ humidified incubator.

Seeding Cells from Frozen Vial:

1. Remove flask and wipe with 70% ethanol. Then transfer to the tissue culture hood. Refer to protocols, flow diagrams and videos for more detail.
2. Transfer the supernatant contents to a sterile centrifuge tubes, this will contain detached cells. Centrifuge at 100g for 7 minutes to obtain cell pellet. Plate cells in a T25 flask pre-coated with ECM.
3. Add pre-warmed Human Bone Marrow Derived Endothelial Stem Cells Growth Media to flask that has attached cells and incubate for two hours in the incubator. After 2 hours wash the cells with 1X PBS, remove 100% 1X PBS solution.
4. Trypsinize the attached cells with 3ml 1X Trypsin EDTA for 2-3 minutes. At the end of Trypsinization add equal volume of complete growth media with serum. Centrifuge the cells at 100g for 7 minutes to obtain cell pellet.
5. Add the 500ul of cells to T25 flask pre-coated with Human Bone Marrow Derived Endothelial Stem Cells Extra-cellular Matrix with 7ml of Human Bone Marrow Derived Endothelial Stem Cells Complete Growth Medium. Depending on the size of the pellet you may utilize a new flask or you may combine the cells to the flask that had the supernatant cells. For T75 Flask add 500ul cell suspension and 10 ml complete growth media.
6. Incubate the cells in the T25or T75 flasks at 37°C in a 5% CO₂ humidified incubator. Perform 100% Media Change every 24 to 48 hours.
7. Medium renewal every other day or 2-3 days, sub-culturing ratio: 1:2 or 1:3 depending on the cell density.
8. Refer to protocols, flow diagrams and videos for more detail. <http://celprogen.com/tech.htm>

Seeding cells from Plated Tissue culture flasks:

For plated cells in tissue culture flask, upon receipt of the cells wipe the flask with 70% ethanol and transfer to sterile tissue culture hood. In the tissue culture hood remove the media from the cells and wash the cells with 1X PBS sterile solution, for 2-3 minutes, remove the 1X PBS solution and then Trypsinize. After Trypsinization



3914 DEL AMO BLVD.
SUITE 901
TORRANCE, CA 90503
www.celprogen.com

Phone: 310 542 8822
Fax: 310 542 8028
Email: info@celprogen.com
stemcells@celprogen.com

of the Cells neutralize the Trypsin with equal volume of Human Bone Marrow Derived Endothelial Stem Cells Complete Growth Media with Serum and collect the Cell suspension in sterile conical centrifuge tube in the tissue culture hood. Centrifuge the cell suspension at 100g for 7 minutes in centrifuge. Plate cells 5×10^5 cells per pre-coated flasks with Human Bone Marrow Derived Endothelial Stem Cells Extra-cellular Matrix for Expansion in Human Bone Marrow Derived Endothelial Stem Cells Complete Growth Medium.

Freezing Medium: Available for purchase Cat# M36096-26FM

Trypsin: Available for purchase Cat# T1509-014

IX PBS: Available for Purchase Cat# P1408-013

Storage temperature: Liquid nitrogen vapor phase

Product Orders: Before submitting an order you will be asked to read and accept the terms and conditions of Celprogen's Material Transfer Agreement (MTA).

Permits/Forms: In addition to the MTA mentioned above, other CELPROGEN and/or regulatory permits may be required for the transfer of this CELPROGEN material. Anyone purchasing CELPROGEN material is ultimately responsible for obtaining the permits.

Notices

& Disclaimers: *CELPROGEN products are intended for laboratory research purposes only. They are not intended for use in humans.* The product, Human Bone Marrow Derived Endothelial Stem Cells, is established and manufactured by CELPROGEN Inc., and is for Research Use Only. This product is not for re-sale or may not be transferred to a third party prior to written request and approval by CELPROGEN Inc.